



CS 30/15 PRIME M

Cod. PR-CS200-040-M

Oil free reciprocating air compressors with soundproof cabinet, 2 cylinders, with a compressed air FAD of 8.12 @ 72.5 PSI pressure. Receiver capacity 40 L, with internal surface treatment to protect against corrosion. Complete with membrane dryer to ensure maximum air purity. This compressor is available also in [145 PSI VERSION](#), visit the appropriate section!



Maximum simultaneous use: 2-3




Dental Chairs: 3

Medical Device CE 0425



Cod.	Air Dryer	Air capacity	FAD Fad @ 72.5 PSI	Motor	Voltage Hz	Receiver	Noise Level	RPM	Duty Cycle	max P.	Size	net Weight
CS 30/15 PRIME M PR-CS200-040-M	Yes	cfm 8.12	cfm 6	HP 2	230V 50Hz	GAL 10.57	db(A) 53	1360	70 %	PSI 116	inch 20.9x29.5x36	Lb 207.3
CS 30/15 PRIME M PR-CS200-040-M-400	Yes	cfm 8.12	cfm 6	HP 2	400V 50Hz	GAL 10.57	db(A) 53	1470	70 %	PSI 116	inch 19.3x28.3x35	Lb 207.3
CS 30/15 PRIME M 220V - 60Hz PR-CS200-040-M-6-UL	Yes	cfm 9.88	cfm 6.53	HP 2	220V 60Hz	GAL 10.57	db(A) 53	1760	70 %	PSI 116	inch 19.3x28.3x35	Lb 207.3
CS 30/15 PRIME M 110V - 60Hz PR-CS200-040-M-110	Yes	cfm 9.88	cfm 6.53	HP 2	110V 60Hz	GAL 10.57	db(A) 53	1760	70 %	PSI 116	inch 19.3x28.3x35	Lb 207.3

Air delivery performance

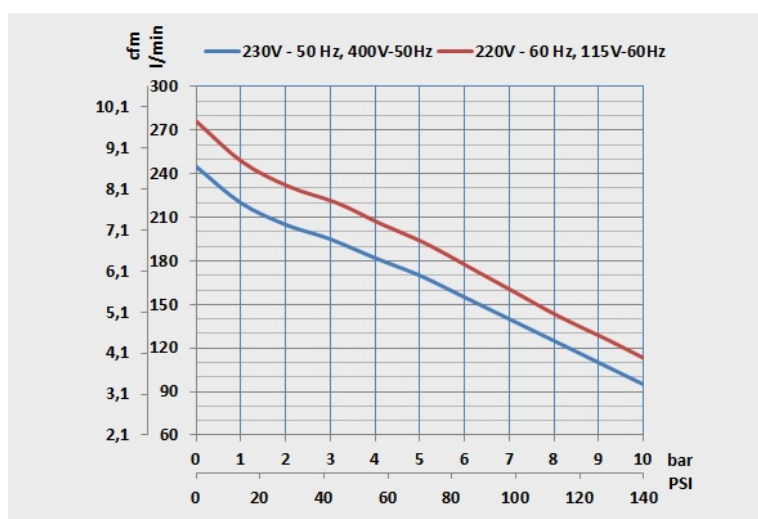

Maintenance kits

Kit.	Description	Type	Hours
TA-FL-026-K-SCP	PRIME 2-3 M filter cartridges kit Cartridge replacement kit for ordinary maintenance of the PURE AIR PRIME 2.3 M range. Comprehends: • 3 suction filters, code: TA-FL-058-C • 1 5 micron filter cartridge, code: TA-FL-001-C-SM • 1 0.01 micron filter cartridge, code: TA-FL-024-C2-SM • 2 o-rings for filters, cod. TA-FL-001-OR • 1 non-return valve pads, code: CO-PV-001 • 1 5 micron filter cartridge (AF30), code: TA-FL-075-C • 1 0.01 micron filter cartridge (AFD30), code: TA-FL-076-C • 2 o-rings for filters, cod. TA-FL-075-OR • 1 AMG150 water separator filter cartridge, code: TA-FL-074-C Replace the filtration cartridges every 2 years.	Ordinary	1.000h/1y
PR-GC-031-K-SCP	PRIME 2 Extraordinary maintenance spare parts kit PRIME 1.5 kW (PRIME 2) versions. comprehends: • 2 suction filters, code: TA-FL-058-C • 4 O-rings, code: MA-PC018 • 2 piston rings , code: MA-MG00001 • 2 piston guide rings, code: MA-MG00003 • 2 expansion rings for the segments, code: MA-MG00024 • 2 complete valve plates, code: MA-MG-02002P-C-P • 1 non-return valve pads, code: CO-PV-001 The extraordinary maintenance intervals recommended for this compressor are 4000 hours/4 years.	Extraordinary	4.000h/4y